

PCR Workstation

TOOLS REQUIRED FOR ASSEMBLY

Electric Hand Drill, Phillips Head Bit, 3/8" Magnetic Hex Driver Bit, Socket Set and Ratchet.

BEFORE ASSEMBLING THE WORKBENCH

Read through the assembly instructions to familiarize yourself with the order in which the parts are assembled.

We suggest that you clear & vacuum the area where the bench is to be assembled.

It is important that the bench is assembled in the same sequence as instructed.



HARDWARE

All assembly hardware is provided.


STEP 1

Flip the desks top upside down on a blanket or other clean surface to avoid scratching the top.

Attach the legs to the modesty panel using (4) 1/4-20 x 1/2" Flange head bolts. The bolts will thread into pre-installed threaded inserts in the sides of the legs




1/4-20 x 1/2"
Flange Head Bolt

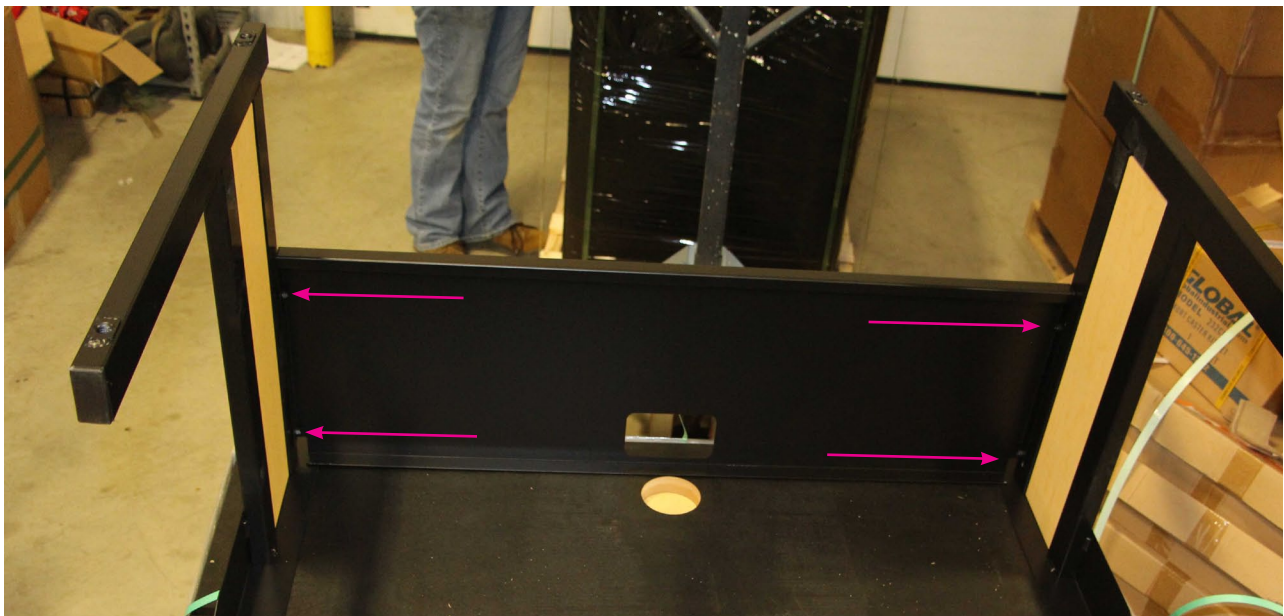


STEP 2

Once legs and modesty panel are attached to each other position the leg mounting holes on the top. The under-side of the top has small prick punches in it. Line the punch up with the holes in the legs. Using (4) #14 x 1" Lag bolts screw into the top by inserting the tip of lag into the pre-punch holes and tightening

Prick Punch

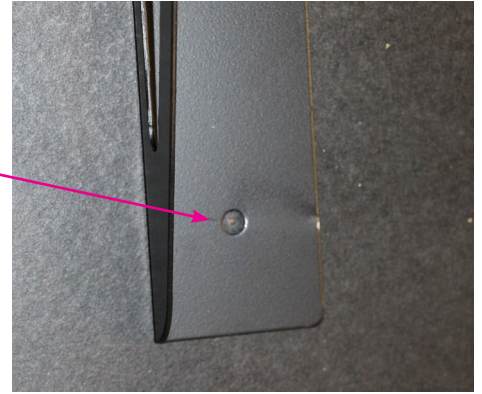

#14 x 1"
Lag Bolt

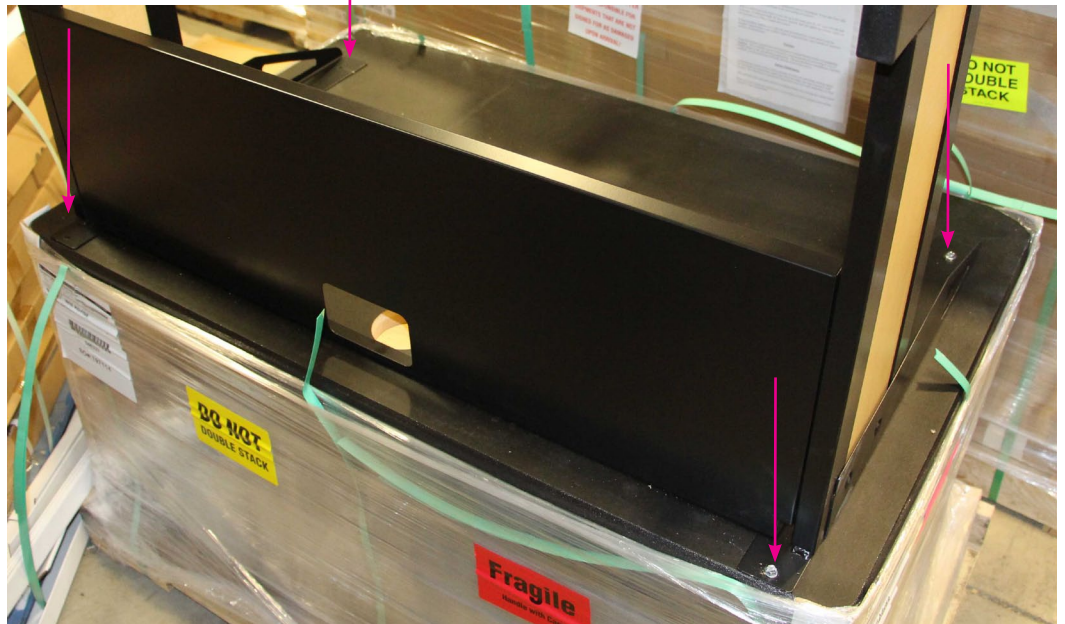

STEP 3

Caster Option:

If ordered thread the single stem casters into the pre-installed nuts on the bottom side of the legs and tighten.

With assistance, flip the desk over right-side-up.

Caster Threaded
Insert

